

stryker®


STRYKER SURGICAL
COMMITTED
TO YOUR SUCCESS.



System 7

Heavy Duty Power Tools

**TOTAL
CONFIDENCE™**

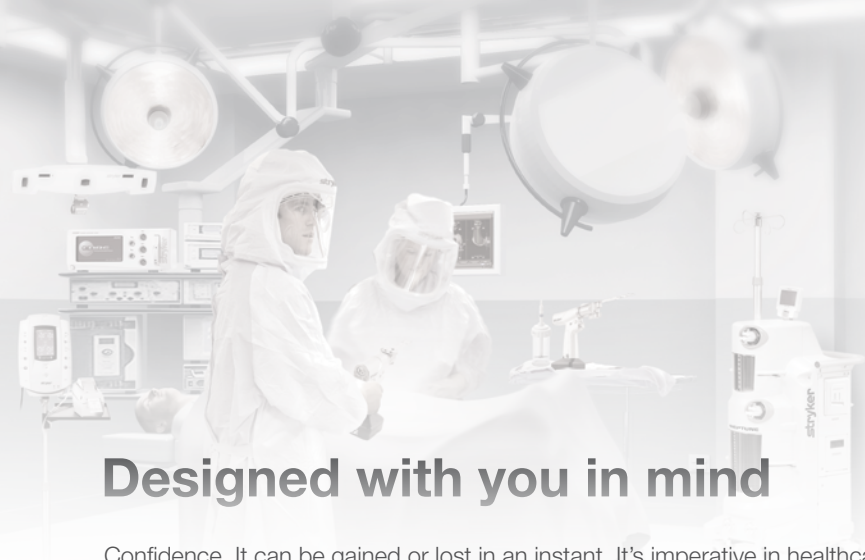

Designed with you in mind

Confidence. It can be gained or lost in an instant. It's imperative in healthcare where decisions can change the course of a life. To make those decisions confidently, you need tools you can depend on. System 7 from Stryker was designed to address the specific needs of your operating room — so you can achieve more, on behalf of your facility and your patients.

Research

System 7 didn't start with us. It started with you: the people working in healthcare facilities, making tough decisions, saving and improving lives. We spent countless hours, and traveled thousands of miles, to discover exactly what you wanted from the next generation of surgical heavy duty power tools — and we listened to every word. The result is an innovative system that meets your demands for quality and reliability.

Development

System 7 has your feedback built in. Every detail has a purpose: to make System 7 the most dependable and highest performing heavy duty power tool system we've ever made.

So you can have the confidence it takes to provide the best care possible.


Stryker System 7.
Total Confidence.

System 7 offers:

- Sagittal Saw is significantly quieter than previous generations
- Stryker SmartLiFE™ Li-Ion Batteries with Exclusive ACP Technology™ mean a longer, more reliable battery life
- Rotary handpiece offers substantially increased torque, providing optimum performance

Sagittal Saw
7208-000-000

Three Decades of Innovation

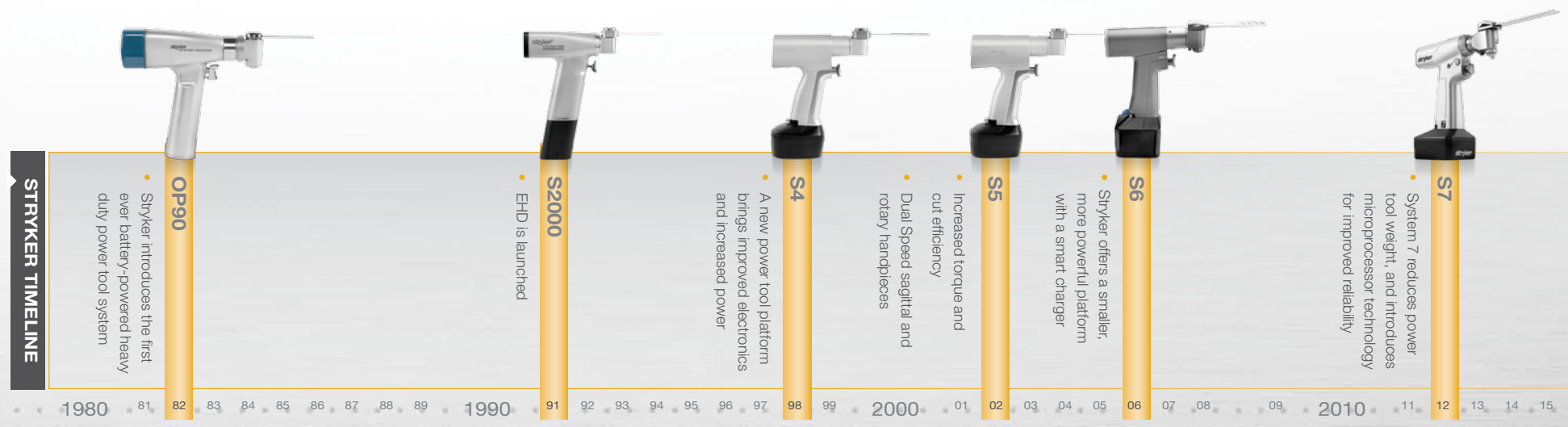

stryker

System 7

TOTAL
CONFIDENCE™

HANDPIECES



Sagittal Saw
7208-000-000

Precision Saw
7209-000-000

Dual Trigger Rotary
7205-000-000

Reciprocating Saw
7206-000-000

Sternum Saw
7207-000-000

Performance

- Improved power-to-weight ratio results in a perfect balance of speed, weight and torque
- Advanced microprocessor continually monitors both the handpiece and battery for optimal performance
- Rotary handpiece offers substantially increased torque, providing optimum performance
- Re-engineered sagittal saw blade mount offers significant noise reduction

Efficiency

- Triathlon attachment eliminates the need for extra attachments
- Unique rotary handpiece design allows all attachments to both drill and ream

Quality & Reliability

- Our most reliable system ever
- All handpieces must pass our 100-Point Quality Check Inspection
- No memory effect enables maximum capacity throughout the life of the battery
- System performs 100 safety checks per second when in use to ensure reliability
- Extended battery shelf life

Ease of Use

- Improved ergonomics:
 - 29% reduction in overall system weight
- Improved sagittal head design offers simple, secure blade loading
- System designed and validated for automated washing
- Rotary handpiece triathlon attachment reduces distance to the surgical site

Standardization

- Backwards compatibility with legacy batteries, attachments and cutting accessories helps reduce inventory and replacement costs



System 7

TOTAL
CONFIDENCE™

When your team members enter the operating room, they need to feel confident that their power tools are ready to meet their needs. That's why System 7 features a newly designed Stryker SmartLiFE Battery with Exclusive ACP Technology. It runs longer, weighs less, offers built-in safety mechanisms and provides proactive warnings – so surgical teams can depend on consistency in every case.



Quality & Reliability

- **Stryker SmartLiFE Batteries with Exclusive ACP Technology**
- No memory effect enables maximum capacity throughout the life of the battery
- Stryker SmartLiFE recyclable batteries provide a green solution
- Laser welded to prevent moisture from entering the battery
- Universal charger identifies overexposure to heat
- Extended battery shelf life

Performance

- Increased run time:
 - Small battery runs 220% longer
 - Large battery runs 183% longer
- Improved ergonomics:
 - 29% reduction in overall system weight

Stryker SmartLiFE™ Batteries with Exclusive ACP Technology™

- Lithium Ion cell technology
- Indicator light proactively alerts OR team when battery is near the end of its charge
- Large and small sizes compatible with all System 7 handpieces

Active Cell Protection (ACP)

Proprietary battery safeguards provide:

- Maximum reliability by regulating voltage and alerting users in the event of overexposure in autoclave
- Unmatched safety by automatically removing voltage when the battery is not on a handpiece or charger



SORN

STERILIZATION

System 7

TOTAL
CONFIDENCE™

Real innovation is often found in small details. It's these kinds of details that make System 7 sterilization cases more cost effective, easier to use and more efficient. So you can achieve positive outcomes for your staff, your facility and your patients.

Efficiency

- Extended battery shelf life
- Bar codes eliminate the need for tracking labels
- Redesigned trays make it easier to load handpieces and attachments in sterilization cases
- Sealed containers minimize SPD processing time

Ease of Use

- Industrial grade sterilization cases are intuitively designed to be durable and easy to use
- No memory effect enables maximum capacity throughout the life of the battery

Cost Control

- Ability to purchase sterilization cases and inserts separately
- Stainless steel sterilization cases are more durable and last longer than previous generations

Automated Washer Compatible

- System designed and validated for automated washing



System 7

Heavy Duty Power Tools


TOTAL
CONFIDENCE™

Ordering Information

System 7 brings you leading-edge technology, backed by a committed support team at Stryker.

Handpieces

- 7208-000-000 Sagittal Saw
- 7209-000-000 Precision Saw
- 7205-000-000 Dual Trigger Rotary
- 7203-000-000 Single Trigger Rotary
- 7206-000-000 Reciprocating Saw
- 7207-000-000 Sternum Saw

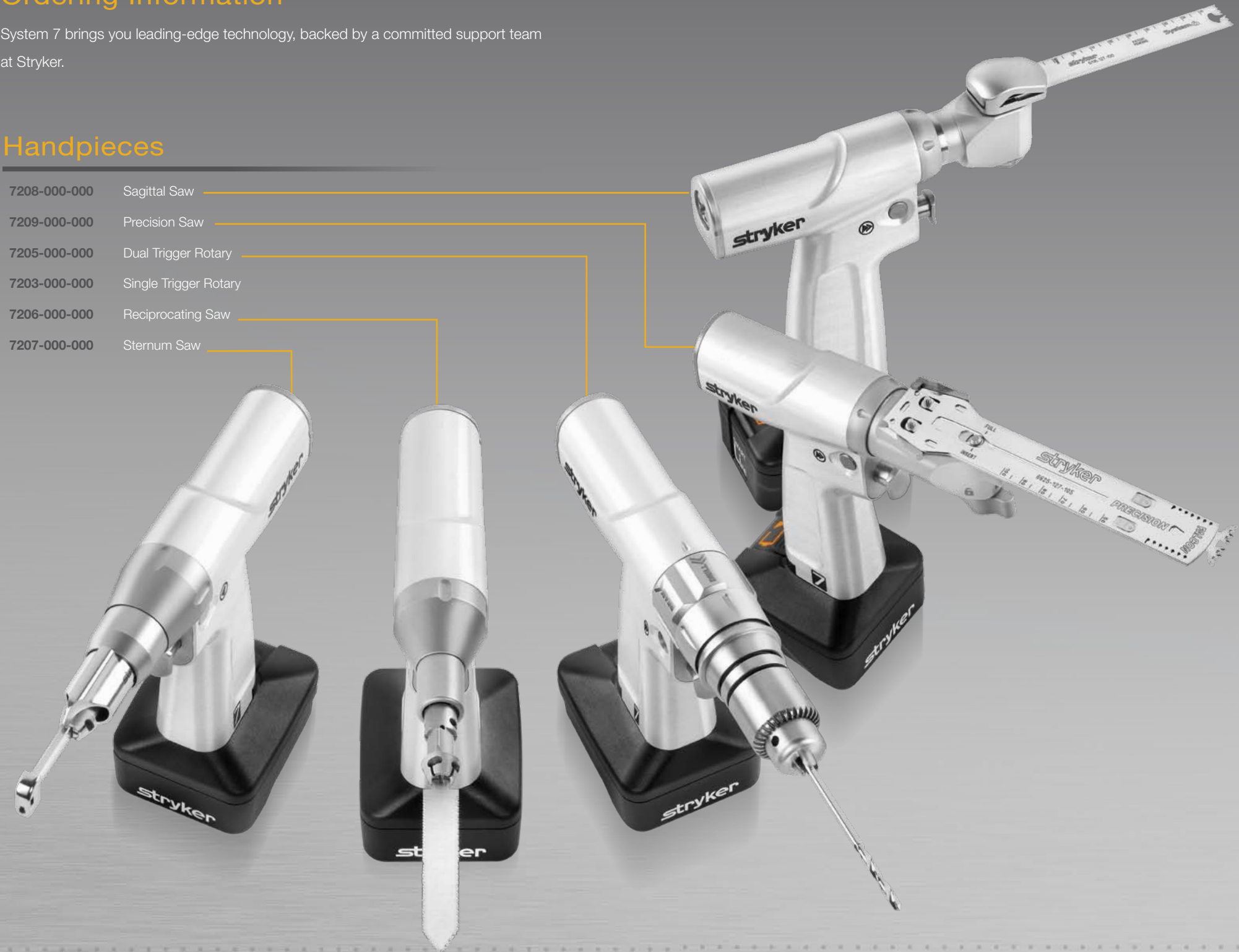

Accessories

- 7203-026-000 Dual Trigger Wire Collet (.7-1.8mm)
- 7203-126-000 Dual Trigger Pin Collet (2.0-3.2mm)
- 7203-226-000 Dual Trigger Pin Collet (3.0-4.2mm)
- 7203-036-000 Single Trigger Wire Collet (.7-1.8mm)
- 7203-136-000 Single Trigger Pin Collet (2.0-3.2mm)
- 7203-236-000 Single Trigger Pin Collet (3.0-4.2mm)
- 6203-110-000 A/O Quick Connect
- 6203-113-000 Hudson
- 6203-131-000 1/4" Keyed Chuck
- 6203-132-000 5/32" Keyed Chuck
- 6203-133-000 1/4" Keyless Chuck
- 6203-134-000 1/8" Keyless Chuck
- 6203-135-000 Hudson Modified Trinkle
- 6203-150-000 Modified Trinkle
- 6203-160-000 Trinkle
- 6203-210-000 Large A/O
- 6203-215-000 DHS / DCS
- 6203-170-000 Triathlon Adapter
- 7207-003-000 Sternum Blade Guard
- 7207-002-000 Sternum Blade Guard Long



Batteries

- 7110-120-000 Universal Battery Charger
- 7212-000-000 Stryker SmartLiFE Small Battery
- 7215-000-000 Stryker SmartLiFE Large Battery



Sterilization Cases

- 7102-450-000 4 HP Sterilization Case Plus Insert Tray
- 7102-454-000 3 HP Sterilization Case Plus Insert Tray
- 7102-453-000 Sternum/Sabo Saw Sterilization Case Plus Insert Tray
- 6102-455-000 Large Aseptic Sterilization Case
- 7102-450-010 4 HP Insert Tray
- 7102-454-010 3 HP Insert Tray
- 7102-453-010 Sternum/Sabo Insert Tray
- 7102-550-000 4 HP Sterilization Container w/Perforated Bottom
- 7102-554-000 3 HP Sterilization Container w/Perforated Bottom
- 7100-553-000 Sternum/Sabo Sterilization Container w/Perforated Bottom





Reconstructive

Hips
Knees
Trauma & Extremities
Joint Preservation
Orthobiologics

Medical & Surgical

Power Tools & Surgical Accessories
Image Guided Navigation
Endoscopy & Arthroscopy
Integrated Communications
Beds, Stretchers & EMS
Sustainability Solutions

Neurotechnology & Spine

Craniomaxillofacial
Interventional Spine
Neurosurgical, Spine & ENT
Neurovascular
Spinal Implants



Stryker Korea
서울시 강남구 삼성동 158-24
동성빌딩 11F
Tel. 02-565-7303
Fax. 02-565-3026
www.stryker.co.kr

본 **stryker**® 카다로그는 전 세계 공용이므로, 본 카다로그에 등재된 제품들 모두가 국내 KFDA 품목허가/신고가 된 것은 아닐 수 있습니다. 이에 제품의 허가 등의 여부를 반드시 사전 확인하시기 바라며, 관련 문의는 **한국스트라이커(주)(대표전화 02.565.7303)**로 연락주시기 바랍니다.

The information presented in this brochure is intended to demonstrate a Stryker product. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.